



## 3-Axes Joystick with Multi-Function Wheel for TANGO Controllers

### Specifications

Number of axes	3 (XYZ)
Interface	HD D-Sub 15-pin male plug connection to TANGO HDI
Material/surface	control panel: membrane keyboard housing: ABS-plastic with haptic coating housing bottom: stainless steel
Dimensions	229 × 139 × 90 mm (L × W × H)
Weight	approx. 1.0 kg


### Order Information

3-Axes Joystick  
 with Multi-Function Wheel  
 Order No.: 00-76-300-0821

### Product Features

- ▶ sensitive, manual positioning
- ▶ side-mounted multi-function wheel allows ergonomic operation
- ▶ various configuration options
- ▶ assignment of the multi-function wheel to any axis
- ▶ free and precise programming of travel speed
- ▶ control of the LED illumination  
LED 100 with multi-function wheel possible
- ▶ 4 freely programmable function keys (e. g. for position capture, velocity select etc.)
- ▶ membrane keyboard with high durability
- ▶ soft touch and insensitive haptic coating on the side surfaces
- ▶ easy cleaning

### Dimensions

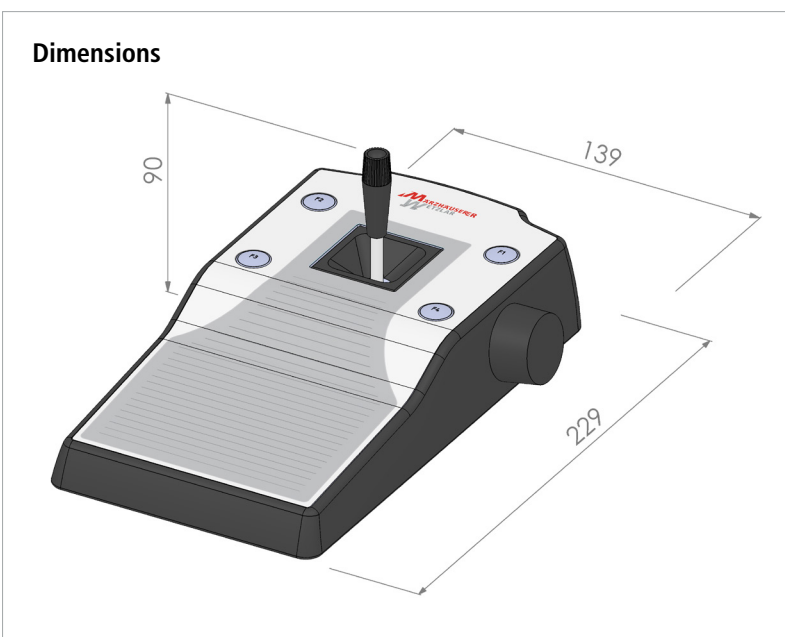

For more accessories for microscopy  
 please visit: [www.marzhauser.com](http://www.marzhauser.com)