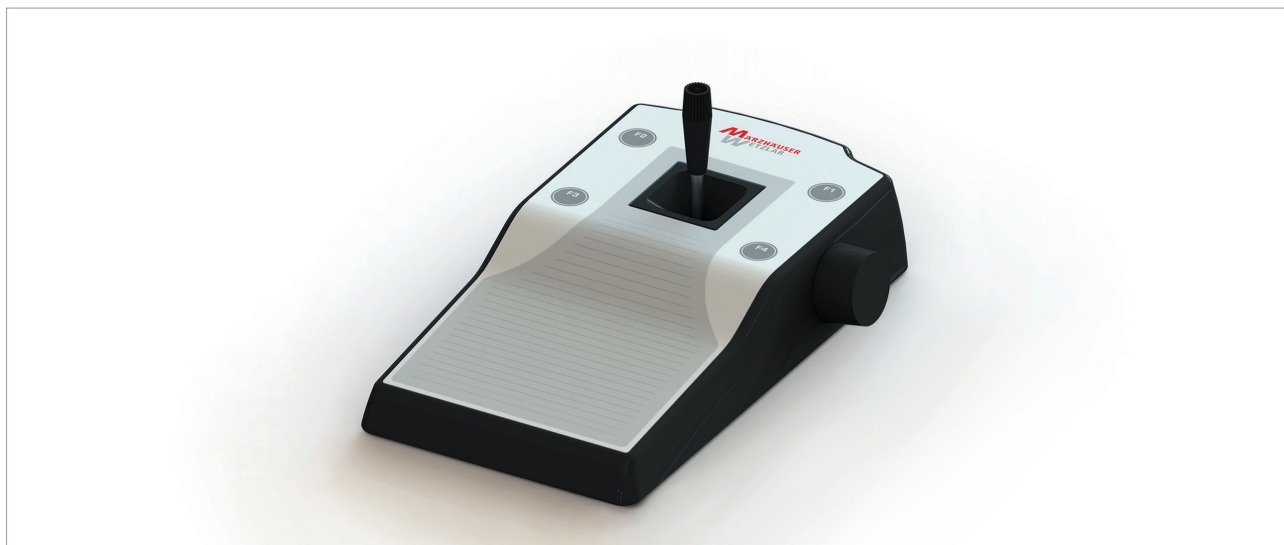




3-Axes Joystick with Multi-Function Wheel for TANGO Controllers

Specifications	analog	digital
Number of axes	3 (XYZ)	
Interface	HD D-Sub 15-pin	Mini USB
Compatibility	TANGO Desktop TANGO PCI-E	TANGO Desktop neo TANGO PCI-E neo TANGO 3 mini ¹ TANGO integrale ¹
Material/surface	control panel: membrane keyboard housing: ABS-plastic with haptic coating housing bottom: stainless steel	
Dimensions	229 × 139 × 90 mm (L × W × H)	
Weight	approx. 1.0 kg	

Order Information

3-Axes Joystick with MF Wheel analog
 Order No.: 00-7604-0310-0000

3-Axes Joystick with MF Wheel digital
 Order No.: 00-7604-1310-0000

Product Features

- ▶ sensitive, manual positioning
- ▶ side-mounted multi-function wheel allows ergonomic operation
- ▶ various configuration options
- ▶ assignment of the multi-function wheel to any axis
- ▶ free and precise programming of travel speed
- ▶ control of the LED illumination
LED 100 with multi-function wheel possible
- ▶ 4 freely programmable function keys (e. g. for position capture, velocity select etc.)
- ▶ membrane keyboard with high durability
- ▶ soft touch and insensitive haptic coating on the side surfaces
- ▶ easy cleaning

For more accessories for microscopy please visit: www.marzhauser.com

Dimensions

